

# REX-KARA II

Complex Beauty Treatment Massager



## REX-KARA II

User Manual & Warranty Certificate




Thanks for purchasing REX-KARA II, a multi-functional massager than enables you to enjoy skin lifting and clearing effects at home. Overall skin care using ultrasonic waves, ion deduction, ion induction, vibration and far infrared ray is available all in one. Please read this manual carefully for long-lasting use and keep it safe with the warranty certificate.

## Safety Instructions

The following instructions are important for your safety. Please read and follow.

	<b>Warning</b>	Circumstances where the user could die or be seriously injured due to misuse.
	<b>Caution</b>	Circumstances where the user could experience disorders or loss.

### Meaning of a Symbol

	Symbol that shows warnings or cautions.	Details about warning and cautions are shown in the pictures or writings near the symbol.
	Symbol that shows restrictions.	Details about restrictions are shown in the pictures or writings near the symbol.
	Symbol that shows compulsions and instructions	Details about compulsions and instructions are shown in the pictures or writings near the symbol.



### Warning

The machine could malfunctions, when used with electric medical devices such as:

- Built-in medical devices, such as an electrocardiograph.
- Life-sustaining devices, such as artificial hearts.

Do not use for patients with:

- Heart Diseases
- Hemophilia
- Malignant Tumors
- Tuberculosis
- Sleeping Disorders, And Is Taking Sleeping Pills
- Blood Pressure Problems
- Contagious Skin Diseases
- Organ Disorders
- Sensation Or Perception Disorders
- Atopic Dermatitis

Do not use for people who are:

- Pregnant
- Under the influence of alcohol
- In menstrual period
- Allergic skin to metal or rubber
- Allergic skin to contacts
- Having inflammation due to cosmetics

Consult your doctor before use if you are:

- Taking hormone drugs
- Having health problems

\* If you notice some troubles on your skin during or after use, consult your doctor before using it again

\* Keep away from children



## Warning


 Unplug	Hold the outlet firmly when plugging in or out. Holding the electric cable instead could start a fire.
	Always unplug when not in use. It could cause fires.
	The outlet is designed for domestic plugs. If it doesn't fit into your plug, do not force it in and contact us.
 No wet hands	Do not touch the outlet or the plug with wet hands. It could cause electric shock.
 Do not disassemble	Do not disassemble, repair or remodel. It could cause malfunction or fire.
 Do not use in water	Like all other electric appliances, do not use in humid conditions. There is danger of electric shock and the machine could fail to work even if the power is off.
	Never use in bathrooms. Keep the device away from water. It could cause mechanical failure, fire, and electric shock.
	Never wash with water. Do not place water or toilet near the machine. It could cause mechanical failure, fire and electric shock.
 Restrictions	Do not damage, bend excessively or cut the electric cord. It could cause fire and electric shock.
	Do not use chargers other than included in the package. It could cause fire and electric shock.
	Do not place inflammable substances near the machine. It could cause explosion or fire.
	Do not use with other electric devices. It could cause machine.
	When during an earthquake or heavy lightning, unplug the machine.
	When the machine has disorders, unplug and stop using.



## Caution


- Do not use longer than the recommended time. The performance will be unsatisfactory.
- Do not keep in humid conditions. It could cause fire and electric shock.
- Do not leave the machine plugged in. It could cause fire.
- When fallen from a high place or damaged, stop using until after inspection. It could cause electric shock and fire.
- When using for a long time, check that the machine is working properly and safely.
- Do not use parts from other machines. It could cause machine failure and electric shock.
- Remove metals from your body when using this machine.
- Always keep your electric neat and tidy.

## Package

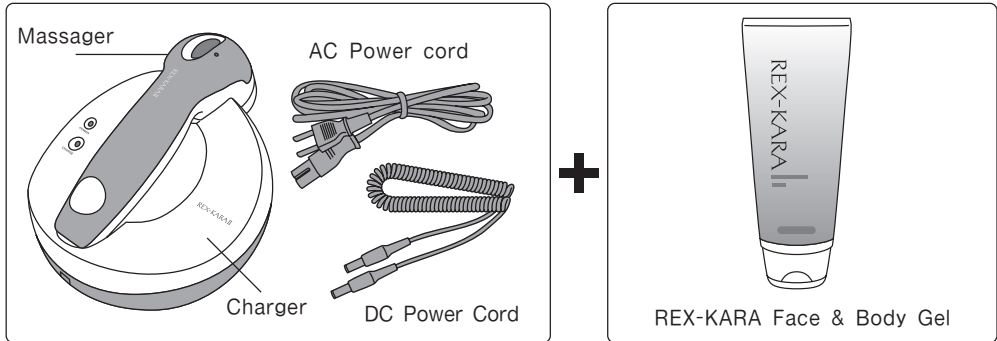


## Name of a Part


### Charger

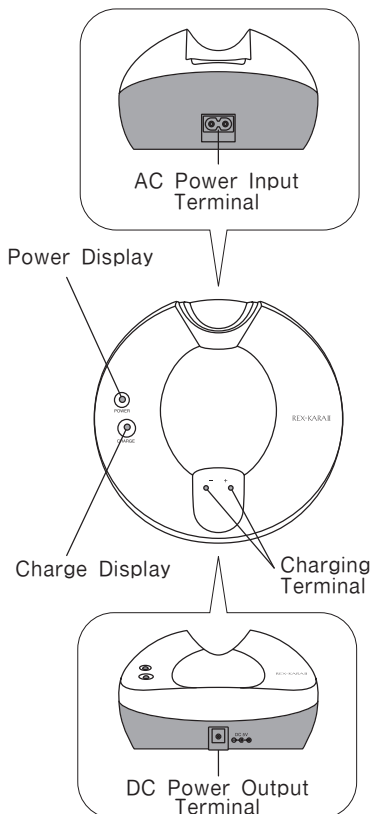

### Handy

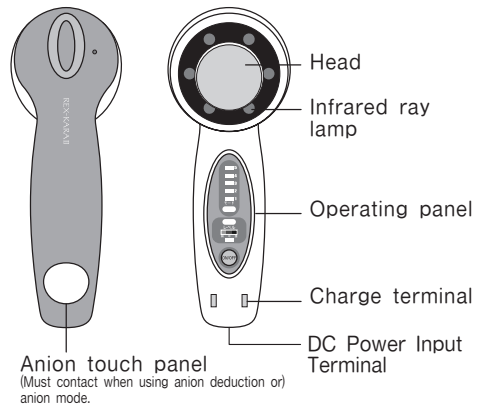




### Acrylic panel

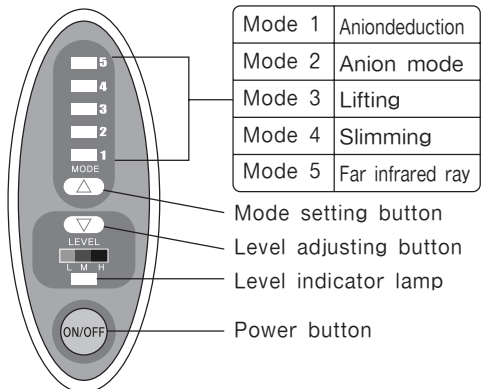

## Before using

### How to charge





1. Plug AC power cord into Ac power input terminal in the charger and plug into outlet at home. (Check that the power display on the charger is red and the charging display is green)
2. Place the massager on the middle of the charger. The charging display turns from green to orange during charging and turns green again when Charging is finished. (Charging time is about 5 hours and using time is about 3 hours).
3. Remove the massager from the charger after charging and use. Also possible to use by connecting DC power cord to the massager and the charger.
4. When the machine can not be charged well, please check if there is waste on the charging terminals both of massager and charger.  
Please wipe with clean cloth to remove waste.

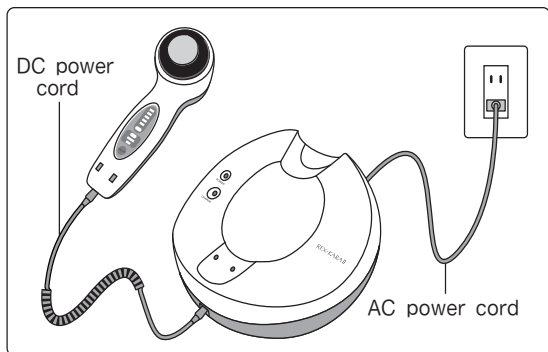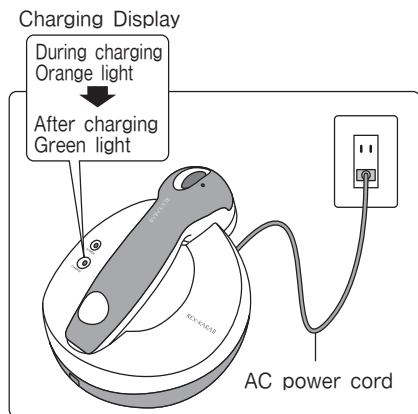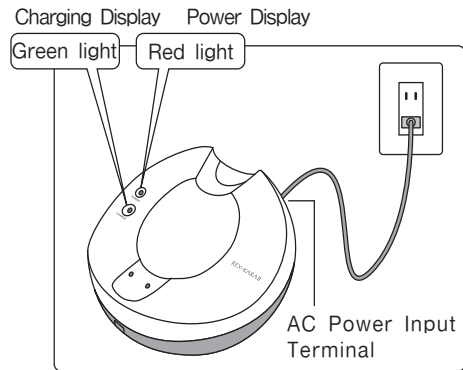

## Before using

This product is a multi-functional massager that enables you to enjoy professional care easily at home. With this one machine, you could experience various skin care using ultrasonic wave, anion deduction, anion mode face care, and far infrared ray. It is easy to use, since all you have to do is select a mode by pressing a button. What's more, you could use it anywhere you want because it is cordless.

### ● Functions of REX-KARA II

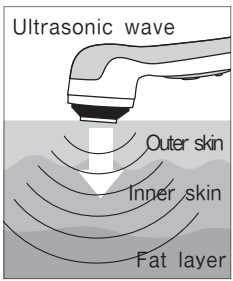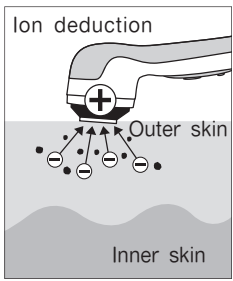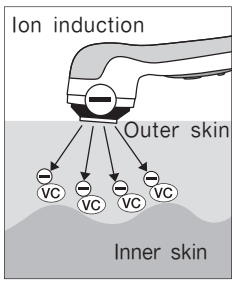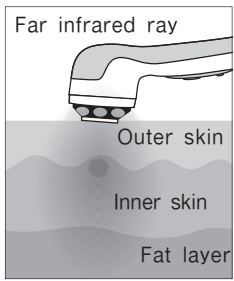 <p><b>Ultrasonic wave</b></p> <p>Outer skin Inner skin Fat layer</p>	 <p><b>Ion deduction</b></p> <p>Outer skin Inner skin</p>	 <p><b>Ion induction</b></p> <p>Outer skin Inner skin</p>	 <p><b>Far infrared ray</b></p> <p>Outer skin Inner skin Fat layer</p>
<p>The wave inaudible to human ear, which has frequency of over 1.7MHz. It produces heat which penetrates deep into inner skin to accelerate metabolism and burns fat</p>	<p>Similar to magnetic attraction, contaminants and wastes inside the skin, which are charged with negative electrons, are pulled out to the positively charged head. Wastes inside the pores are completely removed.</p>	<p>Similar to magnetic repulsion, vitaminC, which is charged with negative electrons, are pushed skin deep using negatively charged head.</p>	<p>Far infrared ray, which penetrates deep inside the body, heats up body from inside. Using this heating action, blood circulation are accelerated and boosts up metabolism.</p>

Titanium coated head

A light, strong and corrosion resistant metal. Used in various areas due to its safety to skin with metal allergies. (Used for jewelry).

## Operating instructions

Before use, wash off cosmetics and contaminants to make your face clean. The effects might be disappointing if there are residues of cosmetics or sebum are left on.

### ● Mode selection

First, turn on the machine by pressing the power button on the operating panel. Next, select a mode from five modes by pressing the mode setting button. (The selected mode lamp will turn green). You will hear an electronic sound whenever a button is pushed. The detailed operating instructions of each mode are available from next page on.

### ● Level selection

After selecting a mode, select level(weak, medium, and strong)by pressing the level adjusting button on the operating panel. The level indicates the power and is divided into weak, medium, and strong depending on the output voltage. Those who want soft care should choose weak and those who seek absolute satisfaction should select strong. Select level that is appropriate for your skin condition. The weak level is indicated by a green lamp, the medium with orange lamp, and the strong level with red lamp. The frequency of ultrasonic wave function, which is 1.7MHz, does not be changed as level.

\* You will hear an electronic sound whenever a button is pushed.

### Mode

Mode 1	Anion deduction
Mode 2	Anion mode
Mode 3	Lifting
Mode 4	Slimming
Mode 5	Far infrared ray

### Level

Weak Green	Medium Orange	Strong Red
---------------	------------------	---------------

## Mode 1 Anion Deduction


Anion Deduction + Ultrasonic Waves

5 Minutes


**Face Care** Completely removes waste inside your pores.



Add plenty of purified water or face lotion to mask sheet or puff and place it on your face

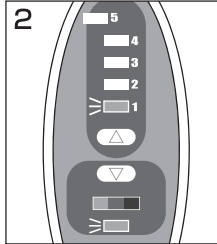

Turn power on and check that the lamp is on after an electronic sound. Mode 1 is selected automatically. Choose level using level adjusting button.

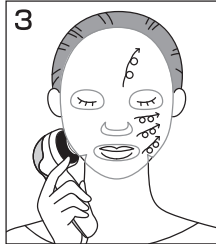

Gently place the handy on top of mask sheet or puff and massage clockwise as shown in the picture. (Massage from center to outside and from bottom to top.)



Remove mask sheet or puff after massaging and wash off remaining waste on your face.

### Cautions when using Mode 1

- \* Always execute mode in contact with the ion touch panel on the back of the handy.
- \* Massage for 20 seconds per spot
- \* Replenish purified water or face lotion to prevent the skin from drying up. (Purified water is sold in pharmacies)
- \* There will be an electronic sound after 5 minutes from starting.

## Mode 2 Anion Induction


Anion Deduction + Ultrasonic Waves + Far-Infrared ray

5 Minutes


**Face Care** A special care that enables treatment substances to penetrate deep inside




Add face lotion to mask sheet or puff and place it on your face

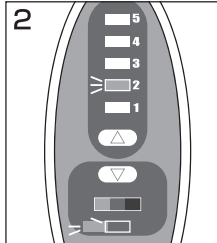

Turn power on and check that the lamp is on after an electronic sound. Select mode 2 by pressing on the mode selecting button one time. Choose level using level adjusting button.

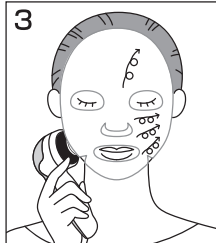

Gently place the massager on top of mask sheet or puff and massage clockwise as shown in the picture. (Massage from center to outside and from bottom to top.)

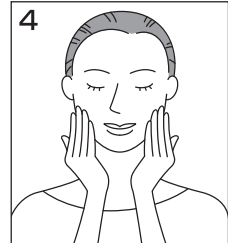

Remove mask sheet or puff after massaging and leave your skin to absorb the remaining substance.

### Cautions when using Mode 2


- \* Always execute mode in contact with the ion touch panel on the back of the massager.
- \* Massage for 20, seconds per spot.
- \* Replenish purified water or face lotion to prevent the skin from drying up.
- \* Far-infrared ray flicker by a regular period.
- \* There will be an electronic sound after 5minutes from starting.
- \* Regular face lotion is also usable.

### Mode **3** Lifting


Ultrasonic + Far - infrared ray

5 Minutes


Face Care Face line lifting effect by ultrasonic waves massage

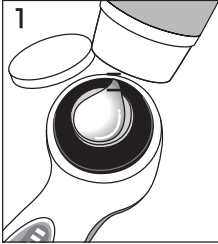

1 Add plenty of gel exclusively used for ultrasonic waves to the head.

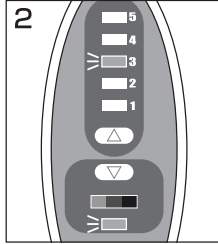

2 Turn power on and check that the lamp is on after an electronic sound. Press mode selecting button twice to select Mode 3. Choose level using level adjusting button

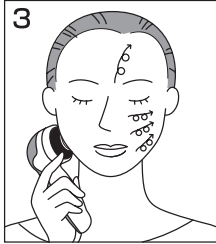

3 Gently place the head on your face and massage clockwise as shown in the picture. (Massage from center to outside and from bottom to top).



4 Remove gel from your face after massaging. (In case of using eye cream instead of gel, give time for your skin to absorb it).

Cautions  
When using  
Mode3


- \* Massage for 20 seconds per spot.
- \* Replenish gel to prevent the skin from drying up.
- \* Far-infrared ray flicker by a regular period.
- \* There will be an electronic sound after 5 minutes from starting.

### Mode **4** Slimming


Ultrasonic Waves

10 Minutes


Body Care Body line slimming effect by ultrasonic waves massage.

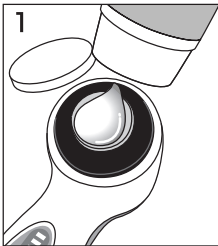


1 Add plenty of gel exclusively used for ultrasonic waves to the head.

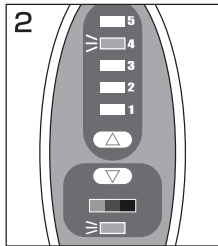

2 Turn power on and check that the lamp is on after an electronic sound. Press mode selecting button three times to select Mode4. Choose level using level adjusting button.

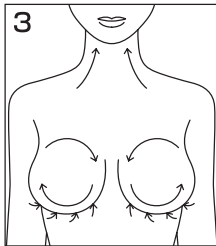

3 For neck, massage from bottom to top and for breasts, massage form bottom to center as shown in the picture.

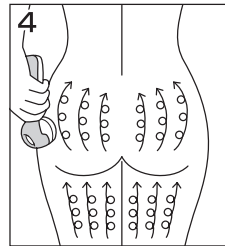

4 For hips, massage from thigh to hips and hips to back in a spiral way.

Cautions  
when using  
Mode4

- \* Massage for 20 seconds per spot.
- \* Replenish gel to prevent the skin from drying up.
- \* There will be an electronic sound after 5 minutes from starting.




Mode **5** Sensitive Care

Far - Infrared ray


10 min

Face Care

Enhances skin condition by generating heat from deep inside your

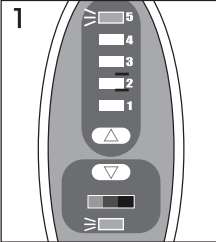

1 Turn power on and check that the lamp is on after an electronic sound. Select mode 5 by pressing on the mode selecting button four times. Choose level using level adjusting button.

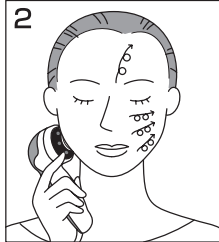

2 Gently place the far-infrared ray lamp on your face and massage clockwise as shown in the picture. (Massage from center to outside and from bottom to top).




3 Wash your face after massaging.




4 After massaging, apply basic skin care agents.

**Cautions**  
When using  
Mode5

- \* Massage for 20 seconds per spot.
- \* Mode 5 could be used on any part of your skin.
- \* There will be an electronic sound after 10 minutes from starting.

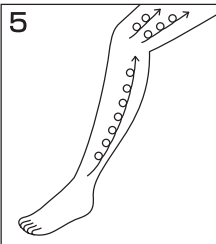

5 For legs, massage from ankle to knee, and knee to thigh in a spiral way.

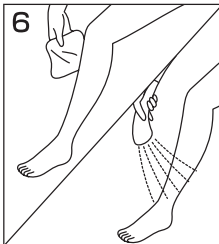

6 Wash off gel from your body after massaging with water.

## Cautions

- Do not use on the eye, eyelid, ear, and pregnant women's stomach.
- Always execute mode 1 and 2 in contact with the ion touch panel on the back of the massager.
- Do not use on dry skin except mode 5.
- Be careful not to liquid or face lotion get in your eyes. If in the case, wash off immediately.
- The massager could become warm during use, but it is not a machine failure.
- Wipe gel off the head after use, to keep it in a clean condition. Do not wash the massager, charger with water.

## FAQ

**Q** I can't feel vibration.

**A** The vibrations of ultrasonic waves are 1.7 million times per second. It is too subtle for you to feel. To check if the massager is working, check the operating sound by putting it close to your ear or check the vibration by applying few drops of water on the head.

**Q** When should I use it?

**A** Use after bath, when your skin is clean and pores are expanded.

**Q** Is it okay to use everyday?

**A** You could REX-KARA II everyday. But the massager is set to work for 5 to 10 minutes, due to dangers of skin damage when over used. Please restrict your time of use to 20 minutes for your face and body respectively. Also, use mode 1, 2, and 3 continuously for outstanding effects. Repeat this total care 2 or 3 times a week.

**Q** Why is it that the mode 1, 2, and 3 are not just one function but a combination of two or more functions?

**A** REX-KARA II has four modes; ultrasonic waves, ion mode, ion deduction, and far infrared rays. We have grouped those that work well together. By selecting a mode, you could use the combined effect without the hassle of choosing several functions.

**Q** Can I use it without recharging?

**A** You could use it without charging, if you connect the charger and the massager with DC Power cord. Use this way when you are unable or have forgotten to charge.

**Q** Is it okay to use gels other than the one included in the package?

**A** Other gels do not guarantee full effect. Please use the gel inside the package when using REX-KARA II.

\*Gels could be ordered separately  
Please make an inquiry to place of purchase.

**Q** Is it okay to wash the massager with water after use?

**A** Do not wash the massager or the charger with water. It could cause failure, fire and electric shock. Wipe off excessive gel with a piece of cloth or tissue.



## Specifications

Material	ABS plastic, Aluminum(titanium, chrome)		
Size	Massager : (W)53 X (L)168 X (H)65mm    AC Power cable : 1.8M Charger : (ø)158 X (H)103mm                DC Power cable : 1.8M		
Weight	Massager : 200g Charger : 300g	Rated frequency	50Hz / 60Hz
Rated Input Charge	AC 100V~240V		
Rated Output Voltage	DC 5.0V/1.5A		
Ultrasonic Wave Frequency	1.7MHz±10% (Vibration : 1.7 million times ±10%)		

## Storage

1. Keep the head of massager clean after use.
2. Do not keep in hot or humid areas.
3. Keep out of children's hands.
4. Unplug the AC power cord when not in use.

## Before inquiring the Customer Relations Center

Please check the below before you inquire for repairs:

1. Is the massager charged? Please place the massager on the charger firmly and recharge.
2. Is the AC Power cord firmly plugged in the outlet?
3. Is the AC Power cord firmly plugged in the AC Power input terminal in the charger?